




Single Use Circular Stapler

Patented Cutting Blade, Breakthrough Innovation in Double Stapling Technology


- Unique staple cutting technology, no sticking of staples after firing
- Safer in low rectal sphincter-preserving procedure
- No damage to anastomosis

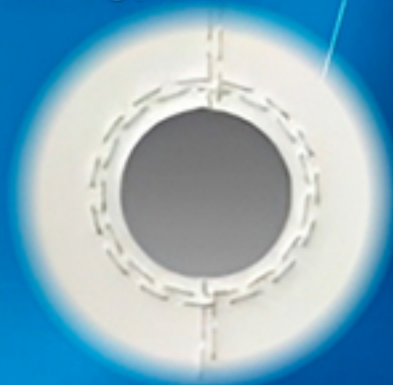

Major Features

- Applicable to end-to-end, end-to-side anastomosis of GI tract and esophagogastric devascularization procedure
- Unique 23mm size, most appropriate for esophageal supra-arch anastomosis
- Precise adjusting of closed staple height, applicable to all tissues
- Gastric tube port in 23 and 26mm models to simplify gastric procedures

Golden Stapler Surgical Co., Ltd.

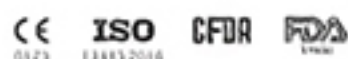

*Advancing Surgical Technology
to Benefit More People*

Patented Staple Cutting Technology

- No sticking of staples after firing
- No damage to anastomosis

5.3mm Open Leg Length

- Staple formation height range: 1.0-2.5mm
- Satisfy stapling of various tissue thickness of gastroenterostomy

Thin and Flat anvil head

- Easy to insert into enteric cavity & remove from anastomosis
- Avoid damage to intestinal wall & anastomosis

Short Distance between Firing and Grabbing Handles

- Simple firing mechanism no exertional effort required

Silicon Coated Handle

- Ergonomic design of great hand touch feeling

Gastric Tube Port

- Load gastric tube in advance, pull it out together with the circular stapler after surgery
- Simplified procedure, reduced duration

5 Different Sizes

- Establish the most appropriate anastomosis
- Allow for various gastroenterostomy

Conjoined Trocar

- Simplify surgical procedures
- Facilitate the lower anterior resection

Plastic Cutting Ring

- Tactile & audible responses
- Indication of safe firing sequence completion

Magnified Green Indication Window with a Gap Setting Scale

- Controllable tissue compression and staple formation height
- Effectively achieve suturing & homeostasis purposes

Reverse Conjoined Trocar

- Applicable for double stapling technique

Code	Size (MM)	Open Staple Height (MM)	Staple Dia. (MM)	Adjustable Closed Staple Height (MM)	Colour	Staple Qty	Circular Knife Dia. (MM)	Staple Rows
CSB21	21	5.3	0.28	1.0-2.5	Jacinth	18	13	2
CSB23	23	5.3	0.28	1.0-2.5	Orange	18	15	2
CSB26	26	5.3	0.28	1.0-2.5	White	20	17	2
CSB29	29	5.3	0.28	1.0-2.5	Blue	24	20	2
CSB32	32	5.3	0.28	1.0-2.5	Green	30	22	2

* Staple from titanium material

Golden Stapler Surgical Co., Ltd.

Add: Building 7A, Jiangsu Wujin Sci-Tech Pioneer Park, 256 Mid Mingxin Road,
Hutang Town, Wujin District, Changzhou, Jiangsu Province, China 213164
Tel: +86 519 85070251 Fax: +86 519 85070252
Web: <http://www.goldenstapler.com>